


BIOCARD™ CELIAC TEST

Sensitivity of Biocard™ Celiac Test was studied with 390 clinical samples and the results were compared to the clinical diagnosis (proven by biopsy).

	Clinical diagnosis	
	+	-
Biocard™	+ 221	6
Celiac Disease	- 5	158

Biocard™ Celiac Test showed 97,4 % sensitivity and 96,9 % specificity.

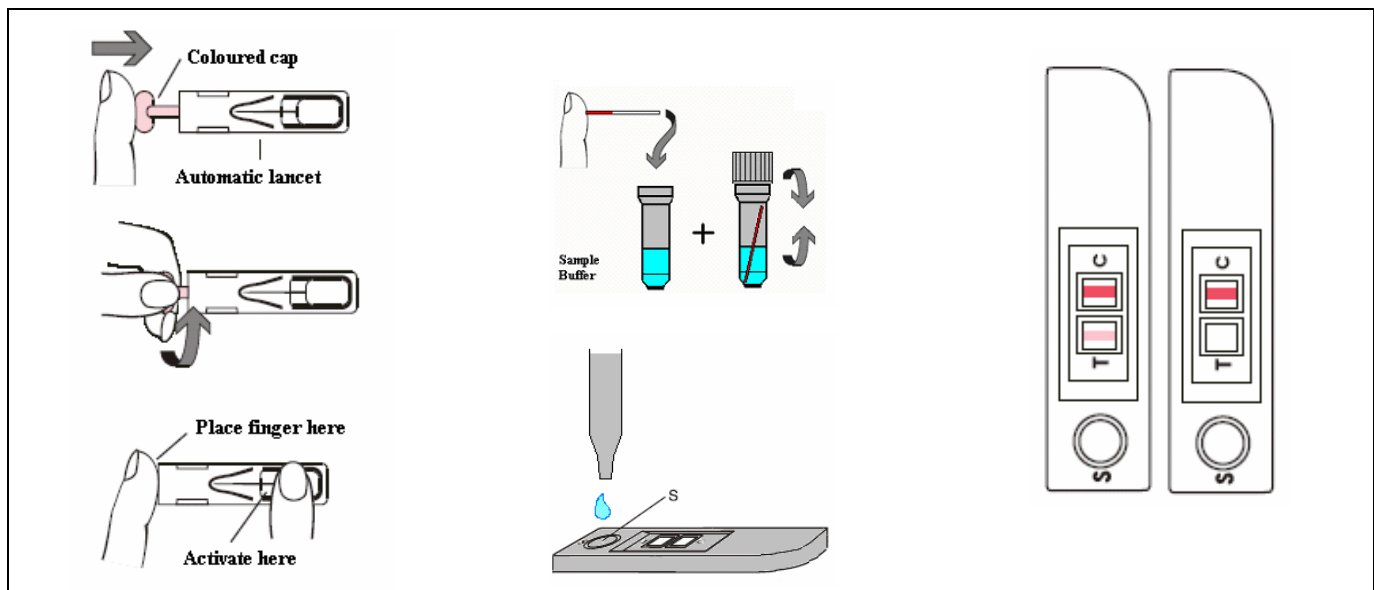

New updated version with total IgA detection. In case of IgA deficiency test will not give control line.

Biocard™ Celiac Test may be used also for monitoring the effects of gluten free diet.

REFERENCES:

1. Aliment Pharmacol Ther 24, 147-154, 2006 T.Raivio et al: Self Transglutaminase-based rapid celiac disease antibody detection by lateral flow method
2. Aliment Pharmacol Ther, 22, 729-737, 2005 I.R. Korponay-Szabo et al: Coeliac disease case finding and diet monitoring by point of care testing
3. Am. J of Gastroenterol, 101, 1-4, 2006 G. Nemeč et al: Looking for celiac disease: Diagnostic accuracy of two rapid commercial assays



Vantaa, Finland
Tel. +358 - 20 155 7510
Fax. +358 - 20 155 7511
www.anibiotech.fi
info@anibiotech.fi